

Necrosis vs Apoptosis Assay Kit

Trial Size (50-100 tests) catalog#KF17371

Standard Size (100-200 tests) catalog#KF17372

FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.

1. INTRODUCTION

Assessment of cellular cytotoxicity levels associated with cytolytic activity of T lymphocytes or natural killer cells is an important aspect of immunobiology related research¹. Early assays developed to assess cytolytic activity utilized the release of ⁵¹Cr from membrane compromised target cells, which were passively loaded with this radioactive indicator prior to exposure to the cytotoxic agent or cells². A major draw-back to these ⁵¹Cr release assays is their reliance on the use of a radioactive isotope with its associated disposal and safe handling issues. Additionally, chromium uptake methods suffer from variabilities in target cell preloading inconsistencies as well the tendency to spontaneously release the chromium isotope into the media in the absence of any cytotoxicity stimulus³.

Assessing potential cytotoxicity properties of chemical and biological agents is a mandatory requirement for the safe distribution of pharmaceuticals, vaccines, or additives associated with food product formulations. Early identification of unintended drug, vaccine, or chemical associated cytotoxicity properties is always an early priority of initial FDA approval testing protocols. With cellular cytotoxicity assessment playing a central role in countless research and environmental safety studies, there is an ever present need for simple, straightforward analysis methods like the Necrosis vs Apoptosis Assay kits developed by Immunochemistry Technologies, LLC (ICT).

Neuromics' Necrosis vs Apoptosis Assay kit simultaneously detects both apoptosis associated cytotoxicity events as well as cell death due to necrosis. Apoptotic cells are identified using ICT's Fluorescent Labeled Inhibitor of Caspases (FLICA) reagent probe⁴⁻⁷. The FAM-FLICA probe covalently binds to active caspase enzymes, which are up-regulated during apoptosis, thus clearly labeling apoptotic cells for subsequent analysis^{8,9}. Non-apoptotic cells will not contain the active caspase enzymes required for FAM-FLICA to remain covalently bound within the cell structure

Loss of the integrity of the cell membrane, indicative of necrosis or late stage apoptosis, is detected using the vital staining dye, 7-aminoactinomycin D (7-AAD), a red fluorescing live/dead stain. This dye easily penetrates cell membrane-compromised cells, binding tightly to GC rich regions of the DNA¹⁰⁻¹⁴. Because 7-AAD alone may not detect cells in the early stages of apoptosis, it is essential to use it in combination with the green-fluorescent FAM-FLICA apoptosis detection reagent¹⁵. Combining these two different types of fluorescent cell-status-indicator reagents within a single test can reveal a significant percentage of cells that are 7-AAD-negative (membrane intact live cells) and yet FAMFLICA positive (apoptotic) (Figures 3-6). FAM-FLICA probes optimally excite between 485-495 nm with maximal emission between 515-525 nm. The vital staining dye 7-AAD can also be efficiently excited at 485-495 nm, but exhibits an optimal emission well into the red fluorescence range (647 nm). This significant difference in fluorescence emission wavelength between the green FAM (carboxyfluorescein) label on the FAM-FLICA probe and the red 7-AAD vital dye simplifies flow cytometer gating and compensation. The FAM-FLICA probe (apoptosis) is monitored on the FL-1 channel, while 7-AAD (necrosis) is monitored on FL-3. Combining the use of ICT's FAMFLICA apoptosis detection probe with a membrane integrity dye like 7-AAD makes it easy to distinguish between necrosis and apoptosis within a single sample

2. KIT CONTENTS

Trial size kits contain:

- 2 vials of FAM-FLICA poly caspase inhibitor reagent #637
- 1 vial of 7-Aminoactinomycin D (7-AAD) vital dye, 0.26 mg vial #6163
- 1 bottle of 10X Apoptosis Wash Buffer (60 mL) #634
- 1 bottle of Fixative (6 mL) #636

FOR RESEARCH USE ONLY

NEUROMICS' REAGENTS ARE FOR IN VITRO AND CERTAIN NON-HUMAN IN VIVO EXPERIMENTAL USE ONLY AND NOT INTENDED FOR USE IN ANY HUMAN CLINICAL INVESTIGATION, DIAGNOSIS, PROGNOSIS, OR TREATMENT. THE ABOVE ANALYSES ARE MERELY TYPICAL GUIDES. THEY ARE NOT TO BE CONSTRUED AS BEING SPECIFICATIONS. ALL OF THE ABOVE INFORMATION IS, TO THE BEST OF OUR KNOWLEDGE, TRUE AND ACCURATE. HOWEVER, SINCE THE CONDITIONS OF USE ARE BEYOND OUR CONTROL, ALL RECOMMENDATIONS OR SUGGESTIONS ARE MADE WITHOUT GUARANTEE, EXPRESS OR IMPLIED, ON OUR PART. WE DISCLAIM ALL LIABILITY IN CONNECTION WITH THE USE OF THE INFORMATION CONTAINED HEREIN OR OTHERWISE, AND ALL SUCH RISKS ARE ASSUMED BY THE USER. WE FURTHER EXPRESSLY DISCLAIM ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. -V2/08/2012

www.neuromics.com

Neuromics Antibodies • 5325 West 74th Street, Suite 8 • Edina, MN 55439
phone 866-350-1500 • fax 612-677-3976 • e-mail: pshuster@neuromics.com

Standard size kits contain:

- 4 vials of FAM-FLICA poly caspase inhibitor reagent #637
- 2 vials of 7-Aminoactinomycin D (vital dye), 0.26 mg vial #6163
- 2 bottles of 10X Apoptosis Wash Buffer (60 mL) #634
- 1 bottle of Fixative (6 mL) #636

3. STORAGE

Store the unopened kit and each unopened component at 2-8°C until the expiration date. Once reconstituted with DMSO, use FAM-FLICA and 7-AAD immediately, or store at ≤-20°C for 6 months protected from light and thawed no more than twice during that time.

4. SAFETY DATA SHEETS (SDS)

Safety data sheets are available upon request

5. RECOMMENDED MATERIALS

- DMSO, up to 1,000 µL (50 µL per vial to reconstitute FAMFLICA, 260 µL per vial to reconstitute 7-AAD, and more to create controls)
- DiH₂O, up to 1,080 mL (540 mL per bottle to dilute 10X Apoptosis Wash Buffer)
- Phosphate buffered saline (PBS) pH 7.4, up to 100 mL, to dilute FAM-FLICA and handle cells
- FBS and/or BSA to add to the buffer when handling cells
- Cultured cells treated with the experimental conditions ready to be labeled
- Reagents to induce caspase activity and create controls for FAMFLICA staining (See Figure 1), such as staurosporine (catalog #6212) or camptothecin (catalog #6210)
- 90% ETOH or 3% formaldehyde to create live/dead controls for 7-AAD staining (See Figure 2)
 - Hemocytometer
 - Centrifuge at < 200 g
 - 15 mL polypropylene centrifuge tubes (1 per sample)
 - Ice bath

6. DETECTION EQUIPMENT

The assay can be analyzed with a:

- Fluorescence microscope
- Flow cytometer Use filter pairings that best approximate these settings:


- FAM-FLICA optimally excites at 488-492 nm and has a peak emission at 515-535 nm (use FL-1 channel).
- 7-AAD optimally excites at 546 nm. It has a peak emission at 647 nm (use FL-3 channel), but displays adequate emission properties when excited within a blue (488-492 nm) light source.

7. EXPERIMENTAL PREPARATION

Staining apoptotic cells with FAM-FLICA and 7-AAD can be completed within a few hours. However, FAMFLICA and 7-AAD are used with living cells, which require periodic maintenance and cultivation several days in advance. In addition, once the proper number of cells has been cultivated, time must be allotted for the experimental treatment or apoptosis induction, which typically requires a 2-6 hour incubation at 37°C based on the cell line and concentration.

Figure 1: Apoptosis Controls Stained with FAM-FLICA

ICT's green poly caspase inhibitor reagent, FAM-FLICA, was used to assess apoptosis in Jurkat cells. Using FAMFLICA, non-apoptotic cells (unstained, left side of each histogram) can easily be distinguished from apoptotic cells exhibiting green fluorescence (right side of each histogram). The forward and side scatter graphs are also shown.

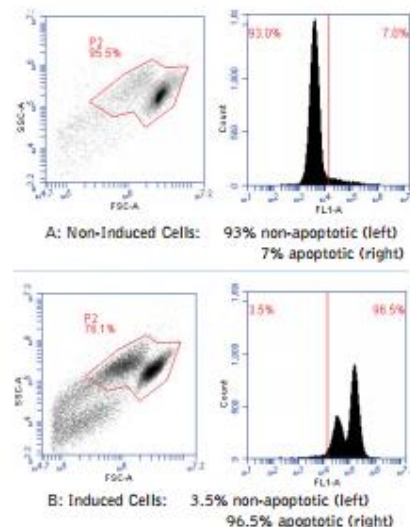

To create a positive control for FAM-FLICA, a population of apoptotic cells is needed (Section 9). If analyzing with a flow cytometer, it will be used to compensate the instrument. In this example, Jurkat cells were grown to 5 x 10⁵ cells/ mL and split into two populations. One population (A, top) was treated with a placebo (DMSO, non-induced) while the other population (B, bottom) was treated with 1 µM staurosporine for 4 hours to induce apoptosis. Cells were stained with FAM-FLICA for 1 hour at 37°C and analyzed using an Accuri C6 flow cytometer in FL-1. Only 7% of non-induced cells (A, right) are apoptotic compared with 96.5% of the induced cells (B, right). (ICT 226:17-19.)

FOR RESEARCH USE ONLY

NEUROMICS' REAGENTS ARE FOR IN VITRO AND CERTAIN NON-HUMAN IN VIVO EXPERIMENTAL USE ONLY AND NOT INTENDED FOR USE IN ANY HUMAN CLINICAL INVESTIGATION, DIAGNOSIS, PROGNOSIS, OR TREATMENT. THE ABOVE ANALYSES ARE MERELY TYPICAL GUIDES. THEY ARE NOT TO BE CONSTRUED AS BEING SPECIFICATIONS. ALL OF THE ABOVE INFORMATION IS, TO THE BEST OF OUR KNOWLEDGE, TRUE AND ACCURATE. HOWEVER, SINCE THE CONDITIONS OF USE ARE BEYOND OUR CONTROL, ALL RECOMMENDATIONS OR SUGGESTIONS ARE MADE WITHOUT GUARANTEE, EXPRESS OR IMPLIED, ON OUR PART. WE DISCLAIM ALL LIABILITY IN CONNECTION WITH THE USE OF THE INFORMATION CONTAINED HEREIN OR OTHERWISE, AND ALL SUCH RISKS ARE ASSUMED BY THE USER. WE FURTHER EXPRESSLY DISCLAIM ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. -V2/08/2012

www.neuromics.com

Neuromics Antibodies • 5325 West 74th Street, Suite 8 • Edina, MN 55439
phone 866-350-1500 • fax 612-677-3976 • e-mail: pshuster@neuromics.com

Create cell populations, such as:

- a. Cells that were exposed to the experimental treatment.
- b. A negative control population of cells that received a placebo treatment.

As FAM-FLICA detects the presence of catalytically active forms of caspase enzymes, plan the experiment so that FAMFLICA will be diluted and administered at the time when caspases are expected to be activated in the cells. Culture cells to a density optimal for the specific experimental conditions or apoptosis induction protocol. Cell density should not exceed 106 cells/mL. Cells cultivated in excess of this concentration may begin to naturally enter apoptosis. An initial experiment may be necessary to determine when and how much FAM-FLICA to use as the resulting positive signal is a direct measurement of caspase activity occurring during the incubation period.

8. CONTROLS

Create experimental samples and control cell populations:

- a. Treated experimental population(s): cells exposed to the experimental condition(s).
- b. Negative control: non-treated cells grown in a normal cell culture environment.
- c. Positive control for FAM-FLICA: cells induced to undergo apoptosis using a known caspase activation protocol (Section 9; Figure 1).
- d. Positive control for 7-AAD: cells treated with ethanol or formaldehyde to create a dead cell population (Section 10 and Figure 2).


The induced positive cell population and negative control cell population tubes should come from a common pool of cells and contain similar quantities of cells. Create negative controls by culturing an equal volume of non-induced cells for every labeling condition. For example, if labeling with FAM-FLICA and 7-AAD stain, make 8 populations:

- 1&2. Unlabeled: induced and non-induced
- 3&4. FAM-FLICA-labeled: induced and non-induced
- 5&6. FAM-FLICA-labeled and 7-AAD-labeled: induced and non-induced
- 7&8. 7-AAD-labeled: induced and non-induced

Figure 2: Live and Dead Controls Stained with 7-AAD

Detection of cell membrane integrity loss, indicative of necrosis or late stage apoptosis, can be done using the red fluorescent live/dead stain, 7-AAD. This vital dye works by penetrating cell

membrane-compromised cells and tightly binding to GC rich regions of the DNA. Using 7-AAD, live cells (unstained, left side of each histogram) can easily be distinguished from dead/membrane compromised cells exhibiting red fluorescence (right side of each histogram). To create a positive control for 7-AAD, a population of dead or killed cell is needed (Section 10). If analyzing with a flow cytometer, it will be used to compensate the instrument. In this example, Jurkat cells were grown to 5×10^5 cells/mL and split into two populations. One population (A) was left untreated while the other population (B) was treated with 90% ethanol for 60 seconds. Cells that were exposed to ethanol (B) were treated with a 5-fold larger volume of PBS to stop the ethanol surface denaturation process. Cells were pelleted by centrifugation (200 x g for 5 minutes) and resuspended in PBS. Cells were then stained with 7-AAD for 10 minutes on ice, and analyzed using an Accuri C6 flow cytometer in FL-3. Only 4% of untreated cells (A) are dead compared with 97.6% of the treated cells (B). Data courtesy of Dr. Kristi Strandberg, ICT 226:30-31.

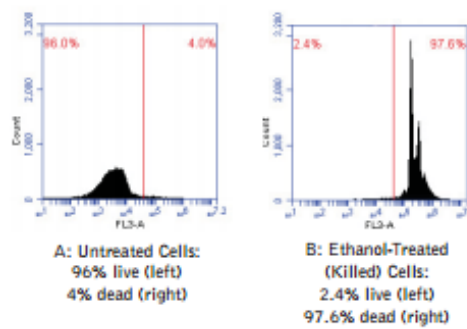

9. APOPTOSIS INDUCTION (FAM-FLICA CONTROL) Prior to commencing the experiment, determine a reproducible method for obtaining a positive control by triggering caspase activity. This process varies significantly with each cell line. For example, apoptosis may be induced with 2-4 $\mu\text{g/mL}$ camptothecin, or 1-2 μM staurosporine for >4 hours (Figure 1).

10. PREPARATION OF DEAD CELLS (7-AAD CONTROL)

Prior to commencing the experiment, determine a reproducible method for obtaining a population of dead/killed cells to use as a positive control for 7-AAD staining. This can easily be achieved using many different techniques. For example, briefly expose cells to 90% ethanol for 30-60 seconds at 37°C, or expose cells to 3% formaldehyde for 30 minutes on ice (Figure 2).

FOR RESEARCH USE ONLY

NEUROMICS' REAGENTS ARE FOR IN VITRO AND CERTAIN NON-HUMAN IN VIVO EXPERIMENTAL USE ONLY AND NOT INTENDED FOR USE IN ANY HUMAN CLINICAL INVESTIGATION, DIAGNOSIS, PROGNOSIS, OR TREATMENT. THE ABOVE ANALYSES ARE MERELY TYPICAL GUIDES. THEY ARE NOT TO BE CONSTRUED AS BEING SPECIFICATIONS. ALL OF THE ABOVE INFORMATION IS, TO THE BEST OF OUR KNOWLEDGE, TRUE AND ACCURATE. HOWEVER, SINCE THE CONDITIONS OF USE ARE BEYOND OUR CONTROL, ALL RECOMMENDATIONS OR SUGGESTIONS ARE MADE WITHOUT GUARANTEE, EXPRESS OR IMPLIED, ON OUR PART. WE DISCLAIM ALL LIABILITY IN CONNECTION WITH THE USE OF THE INFORMATION CONTAINED HEREIN OR OTHERWISE, AND ALL SUCH RISKS ARE ASSUMED BY THE USER. WE FURTHER EXPRESSLY DISCLAIM ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. -V2/08/2012

www.neuromics.com

Neuromics Antibodies • 5325 West 74th Street, Suite 8 • Edina, MN 55439
phone 866-350-1500 • fax 612-677-3976 • e-mail: pshuster@neuromics.com

11. PREPARATION OF FAM-FLICA

FAM-FLICA is supplied as a lyophilized powder that may be slightly visible as an iridescent sheen inside the vial.

Protect from light and use gloves when handling. Because the 30X FAM-FLICA solution must be used immediately, prepare it just before staining.

1. Reconstitute each vial of FAM-FLICA with 50 μ L DMSO to form the 150X stock. The stock solution should be colorless to light yellow. Once reconstituted, it may be stored at $\leq -20^{\circ}\text{C}$ for 6 months protected from light and thawed no more than twice during that time.
 2. Immediately prior to addition to the samples and controls, dilute FAM-FLICA 1:5 by adding 200 μ L PBS to each vial to form the 30X FAM-FLICA solution. Use 30X FAM-FLICA within 30 minutes of dilution into aqueous buffers.
- The recommended volume of 30X FAM-FLICA for flow cytometry is 5 to 10 μ L per 300 μ L of cells at 5×10^5 cells/mL. The recommended volume of 30 X FAM-FLICA for microscopy is 10 μ L per 300 μ L of cells at 5×10^5 cells/mL. These amounts are recommendations, however, the amount needed may vary based on the experimental conditions and the instrument used for analysis. Each investigator should adjust the amount of FAM-FLICA to accommodate the particular cell line and research conditions.

12. PREPARATION OF 1X APOPTOSIS WASH BUFFER

ICT's Apoptosis Wash Buffer (catalog #634, AWB) is an isotonic solution used to wash cells following exposure to FAM-FLICA. It contains mammalian proteins to stabilize cells stained with FAMFLICA, and sodium azide to retard bacterial growth (1X Apoptosis Wash Buffer contains 0.01% w/v sodium azide). Cell culture media containing FBS and other additives may be used to wash cells instead of Apoptosis Wash Buffer.

1. 10X Apoptosis Wash Buffer may form precipitates during cold storage. If this happens, gently warm it until all crystals have dissolved. Do not boil.
 2. Dilute 10X Apoptosis Wash Buffer 1:10 in diH₂O. For example, add 60 mL 10X Apoptosis Wash Buffer to 540 mL diH₂O for a total of 600 mL.
- 1X Apoptosis Wash Buffer may be stored at 2-8 $^{\circ}\text{C}$ and used within 1 week or frozen and used within 6 months.

13. PREPARATION OF 7-AAD

Detection of cell membrane integrity loss, indicative of necrosis or late stage apoptosis, is detected using the red fluorescent live/ dead stain, 7-aminoactinomycin D (7-AAD). This vital dye works by penetrating cell membrane-compromised cells and tightly binding to GC rich regions of the DNA. Because 7-AAD alone may not detect cells in the early stages of apoptosis, it should be combined with the green fluorescent FAM-FLICA reagent to concurrently label apoptotic cells. 7-AAD is supplied as a lyophilized powder that may be slightly visible as a red sheen inside the vial. Protect from light and use gloves when handling.

1. Reconstitute each vial of 7-AAD with 260 μ L DMSO to create a stock concentrate at 1 mg/ mL. Mix by swirling or tilting the vial, allowing the DMSO to travel around the base of the amber vial until the reagent is completely dissolved. At room temperature, the reagent should be dissolved within a few minutes forming a red solution.
2. If storing the stock concentrate for future use, prepare small aliquots (50 μ L) to avoid freeze-thaw cycles. The stock concentrate will be stable for 6 months when protected from light and stored at or below -20 $^{\circ}\text{C}$.
3. When ready to stain cells, use 7-AAD at 1:200 dilution. For example, add 2 μ L per 400 μ L cells.

14. FIXATIVE

ICT's Fixative (catalog #636) is a formaldehyde solution designed to cross-link and aggregate intracellular components. If the stained cell populations cannot be evaluated immediately after labeling with FAM-FLICA, add Fixative at a ratio of 1:10. For example, add 100 μ L Fixative to 900 μ L cells. Never add Fixative until all the staining and final wash steps have been completed. Fixed cells may be stored on ice or at 4 $^{\circ}\text{C}$ for up to 24 hours, protected from light. ICT's Fixative will not interfere with the carboxyfluorescein (FAM) or 7-AAD labels. Do not use absolute ethanol- or methanol-based fixatives as they will inactivate the FAM-FLICA signal. Do not fix cells that will be stained later with 7-AAD (stain the cells before using fixative).

- Danger: Fixative is toxic: danger exists of very serious irreversible effects through inhalation, by contact with skin, or if swallowed. Gloves, protective clothing, and eye wear are strongly recommended. When disposing, flush sink with copious amounts of water; see SDS for further information.

FOR RESEARCH USE ONLY

NEUROMICS' REAGENTS ARE FOR IN VITRO AND CERTAIN NON-HUMAN IN VIVO EXPERIMENTAL USE ONLY AND NOT INTENDED FOR USE IN ANY HUMAN CLINICAL INVESTIGATION, DIAGNOSIS, PROGNOSIS, OR TREATMENT. THE ABOVE ANALYSES ARE MERELY TYPICAL GUIDES. THEY ARE NOT TO BE CONSTRUED AS BEING SPECIFICATIONS. ALL OF THE ABOVE INFORMATION IS, TO THE BEST OF OUR KNOWLEDGE, TRUE AND ACCURATE. HOWEVER, SINCE THE CONDITIONS OF USE ARE BEYOND OUR CONTROL, ALL RECOMMENDATIONS OR SUGGESTIONS ARE MADE WITHOUT GUARANTEE, EXPRESS OR IMPLIED, ON OUR PART. WE DISCLAIM ALL LIABILITY IN CONNECTION WITH THE USE OF THE INFORMATION CONTAINED HEREIN OR OTHERWISE, AND ALL SUCH RISKS ARE ASSUMED BY THE USER. WE FURTHER EXPRESSLY DISCLAIM ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. -V2/08/2012

www.neuromics.com

Neuromics Antibodies • 5325 West 74th Street, Suite 8 • Edina, MN 55439
phone 866-350-1500 • fax 612-677-3976 • e-mail: pshuster@neuromics.com

15. STAINING PROTOCOL FOR SUSPENSION CELLS

Prepare experimental and control cell populations. Ideally, the cell concentration should be 3-5 x 10⁵ cells/mL. The concentration should not exceed 10⁶ cells/mL, as cells cultivated in excess of this concentration may begin to naturally enter apoptosis. Just prior to staining, cells may need to be concentrated to 2-5 x 10⁶ cells/mL as microscopy analysis methods (Section 17) require high cell concentrations. Start with a larger volume of cells at 3-5 x 10⁵ cells/mL (which is a typical density for cell culture) and then concentrate cells and resuspend to 300 µL per sample when ready for FAM-FLICA staining.


1. Expose cells to the experimental or control condition. If analyzing with a flow cytometer, set aside four populations to create instrument control with 7-AAD-positive and 7-AAD-negative cells (Section 10), and FAM-FLICA-induced and FAM-FLICA-non-induced cells (See Figure 5).
2. If analyzing with a fluorescence microscope, concentrate cells to 2-5 x 10⁶ cells/mL just prior to FAM-FLICA staining. Fluorescence microscopy requires an excess of 2 x 10⁶ cells/mL to obtain 5 - 20 cells per image field. Flow cytometry can analyze samples at 3-5 x 10⁵ cells/mL. Sample fluorescent microscopy results are shown in Figure 3.
3. Transfer 290 µL cells into fresh tubes.
4. Add 10 µL 30X FAM-FLICA solution, forming a final volume of 300 µL. If different cell volumes were used, add 30X FAM-FLICA at a ratio of 1:30. Mix by gently flicking the tubes.
 - a. The amount of FAM-FLICA used for flow cytometry applications can be diluted 2-fold for some cell lines and experimental conditions, allowing the end user to test a greater number of samples per kit. If staining a lower concentration, add 5 µL of FAM-FLICA to 295 µL of cells. The amount of FAM-FLICA should be optimized for each cell line and experimental condition.
5. Incubate cells at 37°C protected from light. The incubation period may range from 30 minutes to several hours and should be optimized for each cell line and experimental condition. As cells may settle on the bottom of the tubes, gently resuspend them by swirling cells every 20 minutes to ensure an even distribution of FAM-FLICA throughout the staining process.
6. Centrifuge at 200 x g for 5 minutes, and then discard the supernatant.

7. Add 2 mL 1X Apoptosis Wash Buffer and gently mix.
8. Centrifuge at 200 x g for 5 minutes at RT.
9. Carefully remove and discard supernatants. Gently vortex the pellets to disrupt clumping. Resuspend in 1 mL 1X Apoptosis Wash Buffer and gently mix. Additional wash steps may be required if high FAM-FLICA background signal is observed.

Figure 3: Microscopy Analysis of Jurkat Cells Stained with FAMFLICA and 7-AAD

Jurkat suspension cells were exposed to 1 µM of staurosporine for 4 hours at 37°C to induce apoptosis. Cells were dually stained with the green fluorescent FAM-FLICA poly caspase probe to detect apoptosis via caspase activity, and the red fluorescent vital dye 7-AAD to detect necrosis.

The image shown below in panel A reveals 4 populations of cells: 1) Live, unstained cells, which do not fluoresce. 2) Early stage apoptotic cells fluoresce green with FAM-FLICA. 3) Dually stained green and red fluorescing cells represent the population of Jurkat cells in mid-to-late stage apoptosis; these cells have active caspase enzymes and compromised cell membranes. 4) Necrotic cells fluoresce red. Panel B shows a corresponding differential interference contrast (DIC) image, which reveals cell morphology. In Figure 5, a similar experiment was performed and cells were analyzed via flow cytometry to quantitate each population: 97.1% were apoptotic (Figure 5D, cells in the lower right and upper right quadrants fluoresce green; LR: 46.4% + UR: 50.7%). Microscope images were obtained using an Olympus BH-2 photomicroscope equipped with bright field, DIC, and fluorescence optics. FAM-FLICA and 7-AAD were imaged using a 470-490 nm excitation filter plus >520 nm long pass filter tandem. Data courtesy of Dr. Brian W. Lee (ICT 196:70).

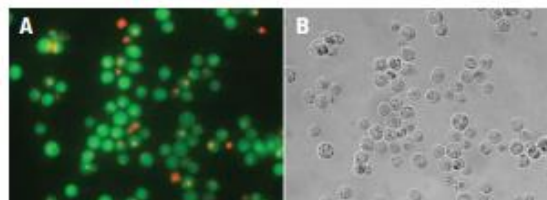

10. Centrifuge cells at 200 x g for 5 minutes at RT.
11. Carefully remove and discard supernatants.
12. Resuspend cells in 400 µL 1X Apoptosis Wash Buffer.
13. Stain with 7-AAD at a final concentration of 5 µg/mL. This can be accomplished by:
 - a. Adding the stock solution directly to the cell culture at 1:200; e.g. add 2 µL stock to 400 µL cell suspension.

FOR RESEARCH USE ONLY

NEUROMICS' REAGENTS ARE FOR IN VITRO AND CERTAIN NON-HUMAN IN VIVO EXPERIMENTAL USE ONLY AND NOT INTENDED FOR USE IN ANY HUMAN CLINICAL INVESTIGATION, DIAGNOSIS, PROGNOSIS, OR TREATMENT. THE ABOVE ANALYSES ARE MERELY TYPICAL GUIDES. THEY ARE NOT TO BE CONSTRUED AS BEING SPECIFICATIONS. ALL OF THE ABOVE INFORMATION IS, TO THE BEST OF OUR KNOWLEDGE, TRUE AND ACCURATE. HOWEVER, SINCE THE CONDITIONS OF USE ARE BEYOND OUR CONTROL, ALL RECOMMENDATIONS OR SUGGESTIONS ARE MADE WITHOUT GUARANTEE, EXPRESS OR IMPLIED, ON OUR PART. WE DISCLAIM ALL LIABILITY IN CONNECTION WITH THE USE OF THE INFORMATION CONTAINED HEREIN OR OTHERWISE, AND ALL SUCH RISKS ARE ASSUMED BY THE USER. WE FURTHER EXPRESSLY DISCLAIM ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. -V2/08/2012

www.neuromics.com

Neuromics Antibodies • 5325 West 74th Street, Suite 8 • Edina, MN 55439
phone 866-350-1500 • fax 612-677-3976 • e-mail: pshuster@neuromics.com

- b. Or by diluting the stock concentrate 1:10 to form the working solution, and then adding the working solution to the cells at 1:20. For example:
 1. Add 50 μ L 7-AAD stock concentrate to 450 μ L PBS or sterile media to form the working solution.
 2. Mix by inverting or vortexing the vial at RT.
 3. Store on ice up to 2 hours.
 4. Add the working solution to the cell suspension at approximately 1:20; e.g. add 25 μ L diluted 7-AAD working solution into 475 μ L cell suspension.
14. Incubate for 10-30 minutes on ice while protected from light.
 - b. To analyze using a microscope, refer to Section 17.
 - a. To analyze by flow cytometer, refer to Section 18.
 - c. Incubate cells with trypsin at 37°C for several minutes. Check frequently to determine when the cells have lifted.
 - d. Neutralize with trypsin inhibitor or cell culture media containing 10-20% FBS.
 - e. Pool cells with any pellets created in Step 2.
3. Centrifuge the collected cells at 200 x g for 5 minutes.
4. Remove all but ~100 μ L supernatant. Resuspend cells in 300-500 μ L cell culture media containing 10-20% FBS. If necessary, count cells and adjust the volume of cell suspension to fit the experiment. Transfer cells into a 15 mL tube.
5. Add 30X FAM-FLICA at 1:30.
 - a. The amount of FAM-FLICA used for flow cytometry applications can be diluted 2-fold for some cell lines and experimental conditions, allowing the end user to test a greater number of samples per kit. If staining a lower concentration, use FAM-FLICA at a dilution between 1:30 and 1:60. The amount of FAM-FLICA should be optimized for each cell line and experimental condition.
6. Incubate 30-60 minutes at 37°C, mixing gently every 10 minutes.
7. Centrifuge at 200 x g for 5 minutes, and then discard the supernatant.
8. Add 2 mL 1X Apoptosis Wash Buffer, and then centrifuge at 200 x g for 5 minutes.
9. Aspirate supernatant and resuspend cells in 2 mL 1X Apoptosis Wash Buffer.
10. Incubate 10 minutes at 37°C to allow any unbound FAMFLICA to diffuse out of the cells.
11. Centrifuge at 200 x g for 5 minutes.
12. Aspirate supernatant and resuspend cells in 400 μ L 1X Apoptosis Wash Buffer.
13. Stain with 7-AAD at a final concentration of 5 μ g/mL. This can be accomplished by:
 - a. Pipetting the stock solution directly into the cell culture at 1:200 (e.g. add 2 μ L 7-AAD stock solution to 400 μ L cell suspension).
 - b. Or by diluting the stock concentrate 1:10 to form the working solution, and then pipetting the working solution into the cells at 1:20. For example:
 1. Add 50 μ L 7-AAD stock concentrate into 450 μ L PBS or sterile media.
 2. Mix by inverting or vortexing the vial at RT.
 3. Store on ice up to 2 hours.

16. STAINING PROTOCOL FOR ADHERENT CELLS

Adherent cells need to be carefully washed to avoid the loss of any cells that round up and come off the plate surface. Loose cells may be harvested from the plate or slide surface and treated as suspension cells, while those remaining adherent to the surface should be washed as adherent cells. If the adherent cells are trypsinized, the loose cells can be recombined with the trypsinized pool; alternatively, the loose cells can be recombined with the adherent portion when the analysis is performed. If growing adherent cells in a plate, the entire plate may be gently spun as part of the wash process to sediment any loose floating cells. If trypsin is necessary, it is preferable to trypsinize prior to FAM-FLICA staining. Avoid trypsinizing cells prior to labeling with a vital dye, like 7-AAD. Cell membranes exposed to trypsin could be transiently permeant to vital dyes for a variable time depending upon the cell line. Cells may be labeled with FAM-FLICA before or after trypsinization (Figure 6).

1. Culture cells in TC-flasks and expose to the experimental conditions.
2. Collect the cells. Note: Apoptotic cells may detach and begin to float into the media. Save and spin to pellet and include these cells in the analysis. If Trypsin is required to lift the adherent cells:
 - a. Rinse the adherent cells with 5-10 mL of PBS to remove any residual trypsin inhibitors present in FBS containing medium.
 - b. Trypsinize by adding enough Trypsin/EDTA to cover the cell surface area.

FOR RESEARCH USE ONLY

NEUROMICS' REAGENTS ARE FOR IN VITRO AND CERTAIN NON-HUMAN IN VIVO EXPERIMENTAL USE ONLY AND NOT INTENDED FOR USE IN ANY HUMAN CLINICAL INVESTIGATION, DIAGNOSIS, PROGNOSIS, OR TREATMENT. THE ABOVE ANALYSES ARE MERELY TYPICAL GUIDES. THEY ARE NOT TO BE CONSTRUED AS BEING SPECIFICATIONS. ALL OF THE ABOVE INFORMATION IS, TO THE BEST OF OUR KNOWLEDGE, TRUE AND ACCURATE. HOWEVER, SINCE THE CONDITIONS OF USE ARE BEYOND OUR CONTROL, ALL RECOMMENDATIONS OR SUGGESTIONS ARE MADE WITHOUT GUARANTEE, EXPRESS OR IMPLIED, ON OUR PART. WE DISCLAIM ALL LIABILITY IN CONNECTION WITH THE USE OF THE INFORMATION CONTAINED HEREIN OR OTHERWISE, AND ALL SUCH RISKS ARE ASSUMED BY THE USER. WE FURTHER EXPRESSLY DISCLAIM ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. -V2/08/2012

4. Add the working solution to the cell suspension at approximately 1:20; e.g. add 25 μ L diluted 7-AAD working solution to 475 μ L cell suspension.
14. Incubate for 10-30 minutes on ice while protected from light. b. To analyze using a microscope, refer to Section 17. a. To analyze by flow cytometer, refer to Section 18.

17. MICROSCOPY ANALYSIS

Follow Section 15 (Steps 1-14), or Section 16 (Steps 1-14).

15. If not viewing immediately, cells may be washed and then fixed for viewing up to 24 hours later.
 - a. Wash sample once with PBS prior to fixing to remove any unbound 7-AAD from the medium.
 - a. Add 30 μ L Fixative. If cells were resuspended in a different volume, add Fixative at a ratio of 1:10.
 - b. Incubate 15 minutes at RT in the dark.
 - c. Place cells onto a microscope slide and allow to dry.
 - d. Briefly rinse cells with PBS.
 - e. Cover cells with mounting media and coverslip. f. Store slides at 2-8°C for up to 24 hours.
16. To view cells immediately, place 1 drop of cell suspension onto a microscope slide and cover with a coverslip.
17. Observe cells under a fluorescence microscope using excitation at 470-490 plus a >520 nm long pass filter. Cells bearing active caspase enzymes covalently coupled to FAM-FLICA appear green. Necrotic or late-stage apoptotic cells containing nucleic acid-bound 7-AAD appear red (Figure 3).

18. FLOW CYTOMETRY ANALYSIS

Follow Section 15 (Steps 1-14), or Section 16 (Steps 1-14).

To address compensation issues and set up the flow cytometer, prepare 4 instrument control populations: Cells stained only with FAM-FLICA (refer to Section 9 and Figure 1):

1. Induced cells.
2. Non-induced cells. Cells stained with only 7-AAD (refer to Section 10 and Figure 2):
3. Live cells.
4. Killed cells (cell membrane compromised).

These controls are needed to adjust the instrument PMT's to separate 7-AAD-positive from 7-AAD-


negative samples and to compensate for bleed-over of the red 7-AAD signal from FL-3 into FL-1. They will also help to clearly differentiate the FAMFLICA-positive population from the FAM-FLICA-negative population and compensate bleed-over of the green FAM-FLICA signal from FL-1 into FL-3.

While setting up the 7-AAD controls (Section 10), continue working with the test samples and experimental controls as prepared in Section 8.

15. Set up the instrument compensation. a. Read the 7-AAD-positive and 7-AAD-negative controls to compensate bleed-over of the red 7-AAD signal from FL-3 into FL-1. b. Read the FAM-FLICA-positive and FAM-FLICA-negative controls to compensate bleed-over of the green FAM-FLICA signal from FL-1 into FL-3.
16. To read the samples for bicolor analysis:
 - a. Measure green FAM-FLICA (FAM = carboxyfluorescein) on the FL-1 channel.
 - b. Measure red fluorescence (7-AAD) on the FL-3 channel.
 - c. Generate a log FL-1 versus log FL-3 dot plot (Figures 5 and 6).
 - d. This will reveal 4 populations of cells: (Figure 4)
 1. Live cells (unstained) only emit low level background fluorescence (lower left).
 2. Cells in early apoptosis fluoresce green with FAMFLICA (lower right).
 3. Cells in late apoptosis are dually stained with FAMFLICA and 7-AAD; they fluoresce green (they have active caspases) and red (the cell membrane is compromised; upper right).
 4. Necrotic membrane-compromised cells fluoresce red (upper left)

Figure 4: Quantitate 4 Cell Populations

1. Live, unstained cells do not fluoresce (LL, lower left).
2. Early stage apoptotic cells fluoresce green with FAM-FLICA (LR, lower right).
3. Dually stained green and red fluorescing cells represent cells in mid to late stage apoptosis; these cells have active caspase enzymes and compromised cell membranes (UR, upper right).
4. Necrotic cells fluoresce red (UL, upper left).

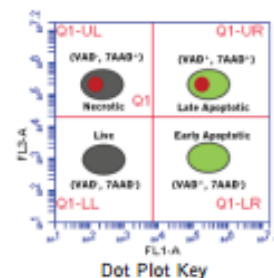

FOR RESEARCH USE ONLY

NEUROMICS' REAGENTS ARE FOR IN VITRO AND CERTAIN NON-HUMAN IN VIVO EXPERIMENTAL USE ONLY AND NOT INTENDED FOR USE IN ANY HUMAN CLINICAL INVESTIGATION, DIAGNOSIS, PROGNOSIS, OR TREATMENT. THE ABOVE ANALYSES ARE MERELY TYPICAL GUIDES. THEY ARE NOT TO BE CONSTRUED AS BEING SPECIFICATIONS. ALL OF THE ABOVE INFORMATION IS, TO THE BEST OF OUR KNOWLEDGE, TRUE AND ACCURATE. HOWEVER, SINCE THE CONDITIONS OF USE ARE BEYOND OUR CONTROL, ALL RECOMMENDATIONS OR SUGGESTIONS ARE MADE WITHOUT GUARANTEE, EXPRESS OR IMPLIED, ON OUR PART. WE DISCLAIM ALL LIABILITY IN CONNECTION WITH THE USE OF THE INFORMATION CONTAINED HEREIN OR OTHERWISE, AND ALL SUCH RISKS ARE ASSUMED BY THE USER. WE FURTHER EXPRESSLY DISCLAIM ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. -V2/08/2012

www.neuromics.com

Neuromics Antibodies • 5325 West 74th Street, Suite 8 • Edina, MN 55439
phone 866-350-1500 • fax 612-677-3976 • e-mail: pshuster@neuromics.com

Figure 5: Flow Cytometry Analysis of Jurkat Suspension Cells to Quantify 4 Populations

Jurkat cells were treated with a placebo (non-induced treatment with DMSO; A and B) or treated with 1 μ M staurosporine for 4 hours to induce apoptosis via caspase activity (C and D). Cells were then dually stained with ICT's FAMFLICA apoptosis reagent and 7-AAD, a vital dye. Cells were analyzed using an Accuri C6 flow cytometer. Compensation was set using cell populations stained individually with either FAM-FLICA or 7-AAD (data not shown). FAM-FLICA was analyzed on FL-1, and 7-AAD was analyzed on FL-3. Forward and side scatter density plots (A, non-induced; C, induced populations) and dual stained density plots (B, noninduced; D, induced populations) are shown. The density plot key is shown in Figure 4.

Flow cytometry can be used to quantitate 4 populations of cells. Live, unstained cells do not fluoresce (lower left, B and D). Early stage apoptotic cells fluoresce green with FAMFLICA (lower right, B and D). Dually stained green and red fluorescing cells represent the population of Jurkat cells in mid to late stage apoptosis; these cells have active caspase enzymes and compromised cell membranes (upper right, B and D). Necrotic cells fluoresce red (upper left, B and D). By including FAM-FLICA in the analysis, cells in early apoptosis (lower right, B and D) can be detected, which 7-AAD cannot detect alone. In the non-induced population (B), only 9.8% of cells were apoptotic (LR: 4.5% + UR: 6.3%) compared with 97.1% of the induced population (D; LR: 46.4% + UR: 50.7%). Cells in these quadrants fluoresce green (ICT 226:17-19). In Figure 3, a similar experiment was performed and images were taken using a fluorescence microscope.


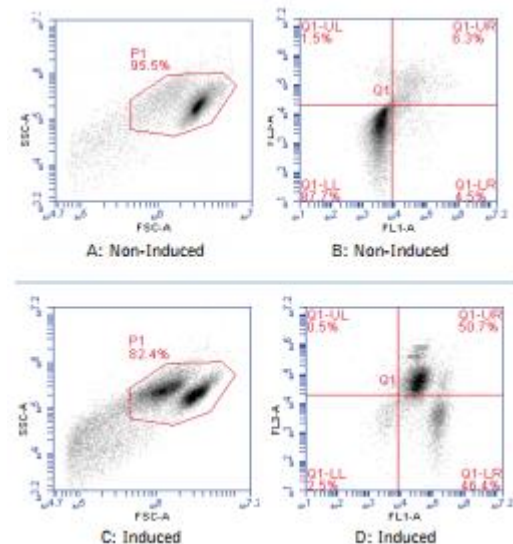


Figure 6: Flow Cytometry Analysis of U-2 OS Adherent Cells to Quantify 4 Populations

U-2 OS cells were treated with a placebo (non-induced, A and B) or treated with 1 μ M staurosporine for 4 hours to induce apoptosis via caspase activity (C and D). Cells were lifted off the tissue culture flask by trypsinization and stained with FAM-FLICA. Cells were then stained with 7-AAD and analyzed by flow cytometry using an Accuri C6 flow cytometer. Compensation was set using cell populations stained individually with either FAM-FLICA or 7-AAD (data not shown). FAM-FLICA was analyzed on FL-1, and 7-AAD was analyzed on FL-3. Forward and side scatter density plots (A, non-induced; C, induced populations) and dual stained density plots (B, non-induced; D, induced populations) are shown. The density plot key is shown in Figure 4. Flow cytometry can be used to quantitate 4 populations of cells. Live, unstained cells do not fluoresce (lower left, B and D). Early stage apoptotic cells fluoresce green with FAM-FLICA (lower right, B and D). Dually stained green and red fluorescing cells represent the population of U-2 OS cells in mid to late stage apoptosis; these cells have active caspase enzymes and compromised cell membranes (upper right, B and D). Necrotic cells fluoresce red (upper left, B and D).

By including FAM-FLICA in the analysis, cells in early apoptosis (lower right, B and D) can be detected, which 7-AAD cannot detect alone. In the non-induced population (B), only 5.9% of cells were apoptotic (LR: 3.6% + UR: 2.3%) compared with 30.5% of the induced population (D; LR: 22.1% + UR: 8.4%). Cells in these quadrants fluoresce green (ICT 226; 17-19).

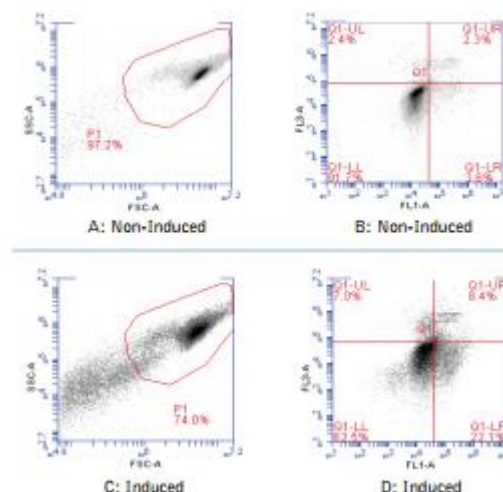

FOR RESEARCH USE ONLY

NEUROMICS' REAGENTS ARE FOR IN VITRO AND CERTAIN NON-HUMAN IN VIVO EXPERIMENTAL USE ONLY AND NOT INTENDED FOR USE IN ANY HUMAN CLINICAL INVESTIGATION, DIAGNOSIS, PROGNOSIS, OR TREATMENT. THE ABOVE ANALYSES ARE MERELY TYPICAL GUIDES. THEY ARE NOT TO BE CONSTRUED AS BEING SPECIFICATIONS. ALL OF THE ABOVE INFORMATION IS, TO THE BEST OF OUR KNOWLEDGE, TRUE AND ACCURATE. HOWEVER, SINCE THE CONDITIONS OF USE ARE BEYOND OUR CONTROL, ALL RECOMMENDATIONS OR SUGGESTIONS ARE MADE WITHOUT GUARANTEE, EXPRESS OR IMPLIED, ON OUR PART. WE DISCLAIM ALL LIABILITY IN CONNECTION WITH THE USE OF THE INFORMATION CONTAINED HEREIN OR OTHERWISE, AND ALL SUCH RISKS ARE ASSUMED BY THE USER. WE FURTHER EXPRESSLY DISCLAIM ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. -V2/08/2012

www.neuromics.com

Neuromics Antibodies • 5325 West 74th Street, Suite 8 • Edina, MN 55439
phone 866-350-1500 • fax 612-677-3976 • e-mail: pshuster@neuromics.com

19. REFERENCES AND CITATIONS

1. Russell, J. H. & Ley, T. J. Lymphocyte-mediated cytotoxicity. *Annu Rev Immunol.* 20, 323-370, Doi:10.1146/annurev.immunol.20.100201.131730 (2002).
2. Brunner, K. T., Mauel, J., Cerottini, J. C. & Chapuis, B. Quantitative assay of the lytic action of immune lymphoid cells on 51-Cr-labelled allogeneic target cells in vitro; inhibition by isoantibody and by drugs. *Immunology.* 14, 181-196 (1968).
3. Jerome, K. R., Sloan, D. D. & Aubert, M. Measurement of CTL-induced cytotoxicity: the caspase 3 assay. *Apoptosis.* 8, 563-571, doi:10.1023/A:1026123223387 (2003).
4. Bedner, E., Smolewski, P., Amstad, P. & Darzynkiewicz, Z. Activation of caspases measured in situ by binding of fluorochrome-labeled inhibitors of caspases (FLICA): correlation with DNA fragmentation. *Exp Cell Res.* 259, 308-313, doi:10.1006/excr.2000.4955 (2000).
5. Ekert, P. G., Silke, J. & Vaux, D. L. Caspase inhibitors. *Cell Death Differ.* 6, 1081-1086, doi:10.1038/sj.cdd.4400594 (1999).
6. Darzynkiewicz, Z., Bedner, E., Smolewski, P., Lee, B. W. & Johnson, G. L. Detection of caspases activation in situ by fluorochrome-labeled inhibitors of caspases (FLICA). *Methods Mol Biol.* 203, 289-299, doi:10.1385/1-59259-179-5:289 (2002).
7. Smolewski, P., Grabarek, J., Halicka, H. D. & Darzynkiewicz, Z. Assay of caspase activation in situ combined with probing plasma membrane integrity to detect three distinct stages of apoptosis. *J Immunol Methods.* 265, 111-121 (2002).
8. Amstad, P. A. et al. Detection of caspase activation in situ by fluorochromelabeled caspase inhibitors. *Biotechniques.* 31, 608-610, 612, 614, passim (2001).
9. Grabarek, J., Amstad, P. & Darzynkiewicz, Z. Use of fluorescently labeled caspase inhibitors as affinity labels to detect activated caspases. *Hum Cell.* 15, 1-12 (2002).
10. Murata, S., Herman, P. & Lakowicz, J. R. Texture analysis of fluorescence lifetime images of nuclear DNA with effect of fluorescence resonance energy transfer. *Cytometry.* 43, 94-100 (2001).
11. Gill, J. E., Jotz, M. M., Young, S. G., Modest, E. J. & Sengupta, S. K. 7-Amino-actinomycin D as a cytochemical probe. I. Spectral properties. *J Histochem Cytochem.* 23, 793-799 (1975).
12. Lecoeur, H., Ledru, E., Prevost, M. C. & Gougeon, M. L. Strategies for phenotyping apoptotic peripheral human lymphocytes comparing ISNT, annexin-V and 7-AAD cytofluorometric staining methods. *J Immunol Methods.* 209, 111-123 (1997).
13. Madhavarao, M. S., Chaykovsky, M. & Sengupta, S. K. N7-Substituted 7-aminoactinomycin D analogues. Synthesis and biological properties. *J Med Chem.* 21, 958-961 (1978).
14. Sengupta, S. K., Tinter, S. K., Lazarus, H., Brown, B. L. & Modest, E. J. 7-substituted actinomycin D analogs. Chemical and growth-inhibitory studies. *J Med Chem.* 18, 1175-1180 (1975).
15. Lee, B. W., Olin, M. R., Johnson, G. L. & Griffin, R. J. In vitro and in vivo apoptosis detection using membrane permeant fluorescent-labeled inhibitors of caspases. *Methods Mol Biol.* 414, 109-135 (2008).

FOR RESEARCH USE ONLY

NEUROMICS' REAGENTS ARE FOR IN VITRO AND CERTAIN NON-HUMAN IN VIVO EXPERIMENTAL USE ONLY AND NOT INTENDED FOR USE IN ANY HUMAN CLINICAL INVESTIGATION, DIAGNOSIS, PROGNOSIS, OR TREATMENT. THE ABOVE ANALYSES ARE MERELY TYPICAL GUIDES. THEY ARE NOT TO BE CONSTRUED AS BEING SPECIFICATIONS. ALL OF THE ABOVE INFORMATION IS, TO THE BEST OF OUR KNOWLEDGE, TRUE AND ACCURATE. HOWEVER, SINCE THE CONDITIONS OF USE ARE BEYOND OUR CONTROL, ALL RECOMMENDATIONS OR SUGGESTIONS ARE MADE WITHOUT GUARANTEE, EXPRESS OR IMPLIED, ON OUR PART. WE DISCLAIM ALL LIABILITY IN CONNECTION WITH THE USE OF THE INFORMATION CONTAINED HEREIN OR OTHERWISE, AND ALL SUCH RISKS ARE ASSUMED BY THE USER. WE FURTHER EXPRESSLY DISCLAIM ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. -V2/08/2012

www.neuromics.com

Neuromics Antibodies • 5325 West 74th Street, Suite 8 • Edina, MN 55439
 phone 866-350-1500 • fax 612-677-3976 • e-mail: pshuster@neuromics.com