


NEUROMICS

ERK2 (MAPK1) Data Sheet

Catalog Number:	RA15001	Host:	Rabbit
Product Type:	Affinity Purified Antibody	Species Reactivity:	Rat, Mouse, Human
Immunogen Sequence:	Recombinant human extracellular signal-regulated kinase-2 (ERK2; accession # NM_138957).	Format:	Liquid 1mg/ml phosphate-buffered saline with 5% trehalose
Applications:	Western blotting: 0.1 ug/ml Immunohistochemistry: 5-15 ug/ml		
Storage:	Dilutions listed only as a recommendation. Optimal dilution should be determined by investigator. The antibody can be stored at 2° - 8° C for 1 month without detectable loss of activity. Antibody can also be aliquotted and stored frozen at -20° C to -70° C in a manual defrost freezer for six months without detectable loss of activity. <i>Avoid repeated freeze-thaw cycles.</i>		

Western blotting

Blotting Buffer	Blocking Solution	Antibody Solution
25 mM Tris, pH 7.4 0.15 M NaCl 0.1% Tween 20	5% nonfat dry milk in Blotting Buffer, Adjust pH to 7.4	5% nonfat dry milk in Blotting Buffer Adjust pH to 7.4

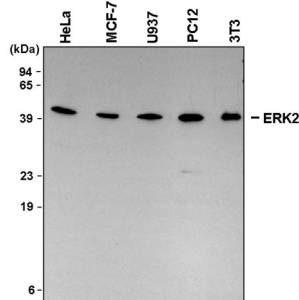

Detection of ERK2 with RA15001.

Lysates from human HeLa and U937 cells and mouse TS1 and NIH/3T3 cells were resolved by SDS-PAGE, transferred to Immobilon membrane and immunoblotted with 0.1 µg/mL goat anti-ERK2,

1. Transfer the electrophoresed proteins to Immobilon membrane (Millipore) and incubate the membrane for 1 hour at room temperature in Blocking Solution.
2. Incubate the membrane overnight at 4° C in Antibody Solution containing 0.1 µg/mL rabbit anti-human/mouse/rat ERK2.
3. Wash the membrane at room temperature for 1 hour with 5 or more changes of Blotting Buffer. Changing the membrane containers often reduces background.
4. Incubate the membrane at room temperature for 1 hour in Antibody Solution containing a HRP-conjugated anti-rabbit IgG secondary antibody.
5. Wash the membrane for 1 hour with 5 or more changes of Blotting Buffer.
6. Detect with ECL Reagent (Amersham).

Cell lysates for Western blottings: To prepare total cell lysates, cells are solubilized in hot 2x SDS gel sample buffer (20 mM dithiothreitol, 6% SDS, 0.25 M Tris, pH 6.8, 10% glycerol, 10 mM NaF and bromophenyl blue) at 2×10^6 - 1×10^7 cells per mL. The extracts are heated in a boiling water bath for 5 minutes and then sonicated with a probe sonicator with 3 - 4 bursts of 5 - 10 seconds each. Samples are diluted with 1x SDS sample buffer to the desired concentration.

Immunohistochemistry

This antibody will detect ERK2 in cells and tissues. The working dilution is 5 - 15 µg/mL.

FOR RESEARCH USE ONLY

NEUROMICS REAGENTS ARE FOR IN VITRO AND CERTAIN NON-HUMAN IN VIVO EXPERIMENTAL USE ONLY AND NOT INTENDED FOR USE IN ANY HUMAN CLINICAL INVESTIGATION, DIAGNOSIS, PROGNOSIS, OR TREATMENT. THE ABOVE ANALYSES ARE MERELY TYPICAL GUIDES. THEY ARE NOT TO BE CONSTRUED AS BEING SPECIFICATIONS. ALL OF THE ABOVE INFORMATION IS, TO THE BEST OF OUR KNOWLEDGE, TRUE AND ACCURATE. HOWEVER, SINCE THE CONDITIONS OF USE ARE BEYOND OUR CONTROL, ALL RECOMMENDATIONS OR SUGGESTIONS ARE MADE WITHOUT GUARANTEE, EXPRESS OR IMPLIED, ON OUR PART. WE DISCLAIM ALL LIABILITY IN CONNECTION WITH THE USE OF THE INFORMATION CONTAINED HEREIN OR OTHERWISE, AND ALL SUCH RISKS ARE ASSUMED BY THE USER. WE FURTHER EXPRESSLY DISCLAIM ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE.

www.neuromics.com

Neuromics Antibodies • 11200 Hampshire Avenue South • Bloomington, MN 55438
phone 507-645-8020 • fax 612-677-3976 • e-mail pshuster@neuromics1.com

03/05v3