



Catalog Number:	MO15108	Host:	Mouse
Product Type:	Protein G Purified, IgG _{2B}	Species Reactivity:	Human, Rat, Mouse
Immunogen Sequence:	Hybridoma resulting from the fusion of a mouse myeloma with B cells obtained from a mouse immunized with purified, <i>E. coli</i> -derived, recombinant human SMC1 (rhSMC1; aa 836 - 1233; Accession # NM_006306).	Format:	Liquid 1mg/ml Solution in phosphate-buffered saline (PBS) with 5% Trehlose
Applications:	Western Blot-0.2-0.5 µg/mL		
Storage:	Dilutions listed as a recommendation. Optimal dilution should be determined by investigator. Antibody can be aliquotted and stored frozen at -20° C to -70° C in a manual defrost freezer for six months without detectable loss of activity. The antibody can be stored at 2° - 8° C for 1 month without detectable loss of activity. <i>Avoid repeated freeze-thaw cycles.</i>		

Application Notes

Specificity

The antibody detects human, mouse and rat SMC1.

Western blot

An antibody concentration of 0.2 - 0.5 µg/mL is recommended.

Protocol

Blotting Buffer

25 mM Tris, pH 7.5
5% nonfat dry milk
0.05% Tween 20

Blocking Solution

5% nonfat dry milk
in Blotting Buffer pH to
7.5

1. Transfer the electrophoresed proteins onto a PVDF membrane and incubate the membrane for 1 hour at room temperature in Blocking Solution.
2. Incubate the membrane for 2 hours at room temperature or overnight at 2 - 8° C in Blocking Solution containing 0.25 - 0.5 µg/mL anti-human/mouse/rat SMC1 antibody.
3. Wash the membrane at room temperature for 30 minutes with 3 or more changes of Blotting Buffer. Changing the membrane containers often reduces background.

FOR RESEARCH USE ONLY

NEUROMICS' REAGENTS ARE FOR IN VITRO AND CERTAIN NON-HUMAN IN VIVO EXPERIMENTAL USE ONLY AND NOT INTENDED FOR USE IN ANY HUMAN CLINICAL INVESTIGATION, DIAGNOSIS, PROGNOSIS, OR TREATMENT. THE ABOVE ANALYSES ARE MERELY TYPICAL GUIDES. THEY ARE NOT TO BE CONSTRUED AS BEING SPECIFICATIONS. ALL OF THE ABOVE INFORMATION IS, TO THE BEST OF OUR KNOWLEDGE, TRUE AND ACCURATE. HOWEVER, SINCE THE CONDITIONS OF USE ARE BEYOND OUR CONTROL, ALL RECOMMENDATIONS OR SUGGESTIONS ARE MADE WITHOUT GUARANTEE, EXPRESS OR IMPLIED, ON OUR PART. WE DISCLAIM ALL LIABILITY IN CONNECTION WITH THE USE OF THE INFORMATION CONTAINED HEREIN OR OTHERWISE, AND ALL SUCH RSKS ARE ASSUMED BY THE USER. WE FURTHER EXPRESSLY DISCLAIM ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE.

03/08v1


www.neuromics.com

Neuromics • 5325 West 74th Street, Suite 8 • Edina, MN 55439
phone 866-350-1500 • fax 612-677-3976 • e-mail pshuster@neuromics1.com

4. Incubate the membrane at room temperature for 1 hour in Blocking Solution containing a 1: 1,000 dilution of goat anti-mouse IgG-HRP
5. Wash the membrane for 30 minutes with 3 or more changes of Blotting Buffer.
6. Detect with WesternGlo Chemiluminescent detection reagents

Cell lysates for Western blotting: To prepare total cell lysates, solubilize cells in 2X SDS gel sample buffer (20 mM dithiothreitol, 6% SDS, 0.25 M Tris, pH 6.8, 10% glycerol, and bromophenyl blue) and sonicate with a probe sonicator using 3 - 4 bursts of 5 - 10 seconds each. Heat extracts in a boiling water bath for 5 minutes and load onto polyacrylamide gels. Samples may be diluted with 1X SDS sample buffer to the desired concentration.

Image: Detection of SMC1 with MO15108. Extracts from exponentially growing U2-OS (human), Balb/c-3T3 (mouse), and normal rat kidney (NRK) cells were prepared, resolved by SDS-PAGE, and transferred to a PVDF membrane. The membrane was immunoblotted with 0.25 µg/mL monoclonal anti-SMC1 antibody.

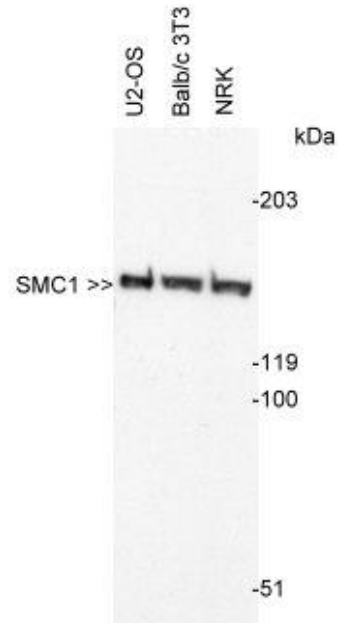

FOR RESEARCH USE ONLY

NEUROMICS' REAGENTS ARE FOR IN VITRO AND CERTAIN NON-HUMAN IN VIVO EXPERIMENTAL USE ONLY AND NOT INTENDED FOR USE IN ANY HUMAN CLINICAL INVESTIGATION, DIAGNOSIS, PROGNOSIS, OR TREATMENT. THE ABOVE ANALYSES ARE MERELY TYPICAL GUIDES. THEY ARE NOT TO BE CONSTRUED AS BEING SPECIFICATIONS. ALL OF THE ABOVE INFORMATION IS, TO THE BEST OF OUR KNOWLEDGE, TRUE AND ACCURATE. HOWEVER, SINCE THE CONDITIONS OF USE ARE BEYOND OUR CONTROL, ALL RECOMMENDATIONS OR SUGGESTIONS ARE MADE WITHOUT GUARANTEE, EXPRESS OR IMPLIED, ON OUR PART. WE DISCLAIM ALL LIABILITY IN CONNECTION WITH THE USE OF THE INFORMATION CONTAINED HEREIN OR OTHERWISE, AND ALL SUCH RISKS ARE ASSUMED BY THE USER. WE FURTHER EXPRESSLY DISCLAIM ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE.

03/08v1

www.neuromics.com

Neuromics • 5325 West 74th Street, Suite 8 • Edina, MN 55439
phone 866-350-1500 • fax 612-677-3976 • e-mail pshuster@neuromics1.com