



Catalog Number: GP10108

Host: Guinea Pig

Product Type: Whole serum

Species: Rat

Immunogen Sequence: VEKQSTDSGAYSIGH-Corresponding to residues 383-397 of the carboxy-terminus of rat P2X3

Reactivity:
Format: Whole Serum (with 0.05% sodium azide). Sent in liquid form

Applications: Immunohistochemistry 1:100 – 1:1,000
Dilutions listed as a recommendation. Optimal dilution should be determined by investigator.

Storage: Store frozen. Aliquot as undiluted serum and immediately place at -20°C. Serum may have become trapped in top of vial during shipping. Centrifugation of vial is recommended before opening. Stable for at least 6 months at -20°C. Repeated freeze/thaw cycles compromise the integrity of the antiserum.

Reference: [Jennifer C. Peleshok, Alfredo Ribeiro-da-Silva. Delayed reinnervation by nonpeptidergic nociceptive afferents of the glabrous skin of the rat hindpaw in a neuropathic pain model.](#) The Journal of Comparative Neurology. Volume 519, Issue 1, pages 49–63, 1 January 2011.

[C.K. Park, J.H. Bae, H.Y. Kim, H.J. Jo, Y.H. Kim, S.J. Jung, J.S. Kim, and S.B. Oh. Substance P Sensitizes P2X3 in Nociceptive Trigeminal Neurons.](#) Journal of Dental Research, October 2010; 89: 1154 – 1159

[James P. Lund, Somayeh Sadeghi, Tuija Athanassiadis, Nadia Caram Salas, François Auclair, Benoît Thivierge, Isabel Arsenault, Pierre Rompré, Karl-Gunnar Westberg, and Arlette Kolta. Assessment of the Potential Role of Muscle Spindle Mechanoreceptor Afferents in Chronic Muscle Pain in the Rat Masseter Muscle.](#) PLoS One. 2010; 5(6): e11131. Published online 2010 June 15. doi: 10.1371/journal.pone.0011131.

Application Notes

Immunohistochemistry:

Antiserum was used on perfusion fixed tissue. Perfusion: 1) calcium-free Tyrode's solution, 2) paraformaldehyde-picric acid fixative, and 3) 10% sucrose in PBS as a cryo-protectant. Desired tissues were dissected and stored overnight in 10% sucrose in PBS.

Slide-mounted tissue sections were processed for indirect immunofluorescence. Slides were incubated with blocking buffer for 1 hour at room temperature. Primary antiserum was diluted with blocking buffer to the appropriate working concentration. Blocking buffer was removed and slides were incubated for 18-24 hours at 4°C with primary antiserum. Slides were rinsed 3 times and then incubated with secondary antibodies for 1 hour at room temperature. Slides were again rinsed 3 times and coverslipped. Staining was examined using fluorescence microscopy.

FOR RESEARCH USE ONLY

NEUROMICS' REAGENTS ARE FOR IN VITRO AND CERTAIN NON-HUMAN IN VIVO EXPERIMENTAL USE ONLY AND NOT INTENDED FOR USE IN ANY HUMAN CLINICAL INVESTIGATION, DIAGNOSIS, PROGNOSIS, OR TREATMENT. THE ABOVE ANALYSES ARE MERELY TYPICAL GUIDES. THEY ARE NOT TO BE CONSTRUED AS BEING SPECIFICATIONS. ALL OF THE ABOVE INFORMATION IS, TO THE BEST OF OUR KNOWLEDGE, TRUE AND ACCURATE. HOWEVER, SINCE THE CONDITIONS OF USE ARE BEYOND OUR CONTROL, ALL RECOMMENDATIONS OR SUGGESTIONS ARE MADE WITHOUT GUARANTEE, EXPRESS OR IMPLIED, ON OUR PART. WE DISCLAIM ALL LIABILITY IN CONNECTION WITH THE USE OF THE INFORMATION CONTAINED HEREIN OR OTHERWISE, AND ALL SUCH RISKS ARE ASSUMED BY THE USER. WE FURTHER EXPRESSLY DISCLAIM ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE.

6/04v1

www.neuromics.com

Neuromics Antibodies • 5325 West 74th Street, Suite 8 • Edina, MN 55439
phone 866-350-1500 • fax 612-677-3976 • e-mail pshuster@neuromics.com

Image: P2X3 staining in rat DRG.


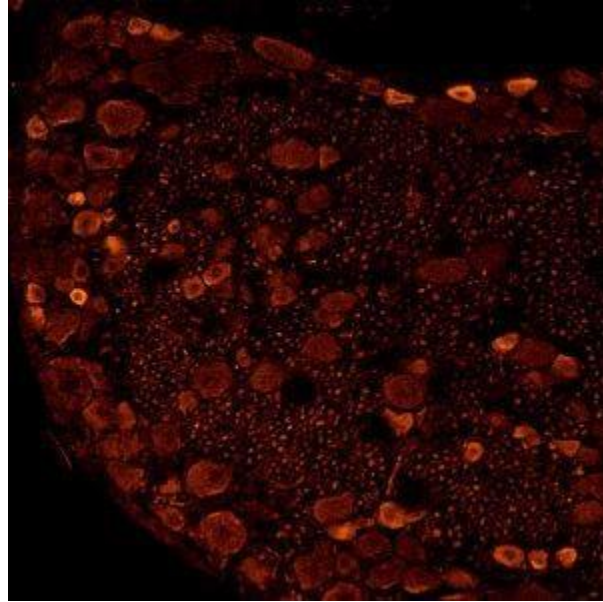


Image: P2X3 staining in Dorsal Horn.

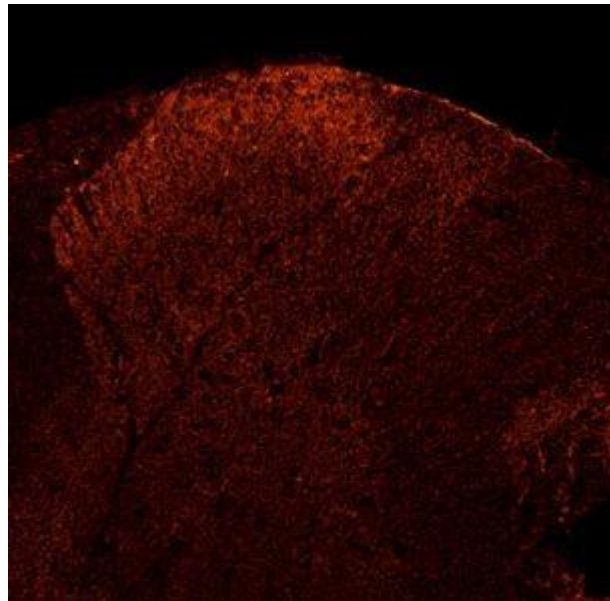

FOR RESEARCH USE ONLY

NEUROMICS' REAGENTS ARE FOR IN VITRO AND CERTAIN NON-HUMAN IN VIVO EXPERIMENTAL USE ONLY AND NOT INTENDED FOR USE IN ANY HUMAN CLINICAL INVESTIGATION, DIAGNOSIS, PROGNOSIS, OR TREATMENT. THE ABOVE ANALYSES ARE MERELY TYPICAL GUIDES. THEY ARE NOT TO BE CONSTRUED AS BEING SPECIFICATIONS. ALL OF THE ABOVE INFORMATION IS, TO THE BEST OF OUR KNOWLEDGE, TRUE AND ACCURATE. HOWEVER, SINCE THE CONDITIONS OF USE ARE BEYOND OUR CONTROL, ALL RECOMMENDATIONS OR SUGGESTIONS ARE MADE WITHOUT GUARANTEE, EXPRESS OR IMPLIED, ON OUR PART. WE DISCLAIM ALL LIABILITY IN CONNECTION WITH THE USE OF THE INFORMATION CONTAINED HEREIN OR OTHERWISE, AND ALL SUCH RSKS ARE ASSUMED BY THE USER. WE FURTHER EXPRESSLY DISCLAIM ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE.

6/04v1

www.neuromics.com

Neuromics Antibodies • 5325 West 74th Street, Suite 8 • Edina, MN 55439
phone 866-350-1500 • fax 612-677-3976 • e-mail pshuster@neuromics.com