# NEUROMICS /



# Neurofilament NF-H

## **Data Sheet**

Catalog Number: MO22187 Host: Mouse

**Species** Human, Rat, Mouse, Cow, Pig, and Mouse Monoclonal IgG **Product Type:** 

Reactivity:

Native NF-H purified from bovine spinal Purified antibody at 1mg/mL in 50% **Immunogen** Format:

PBS, 50% glycerol plus 5mM NaN<sub>3</sub> Sequence:

Immunofluorescence: 1:1,000 Applications:

Immunohistochemistry: 1:1,000 Immunocytochemistry: 1:1,000 Western Blot: 1:10,000

Dilutions listed as a recommendation. Optimal dilution should be determined by investigator.

Storage: Antibody can also be aliquotted and stored frozen at -20° C in a manual defrost freezer for six

months without detectable loss of activity. The antibody can be stored at 2° - 8° C for 1 month


without detectable loss of activity. Avoid repeated freeze-thaw cycles.

### **Application Notes**

#### Description/Data:

Neurofilaments are the 10nm or intermediate filament proteins found specifically in neurons, and are composed predominantly of three major proteins called NF-L, NF-M and NF-H, though other proteins may also be present. NF-H is the neurofilament high or heavy molecular weight polypeptide and runs on SDS-PAGE gels at 200-220 kDa, with some variability across species boundaries. Antibodies to NF-H are useful for identifying axonal processes in tissue sections and in culture. NF-H antibodies can also be useful in visualizing neurofilament accumulations seen in many neurological disorders, such as Amyotrophic Lateral Sclerosis (also known as Lou Gehrig's disease), Alzheimer's disease and following traumatic injury.

Image: Immunohistological analysis of rat cerebellum section stained with mouse mAb to pNF-H, MO22187, dilution 1:2,000 in green, and costained with rabbit pAb to FOX3/NeuN, dilution 1.5,000 in red. Following transcardial perfusion with 4% paraformaldehyde, brain was post fixed for 24 hours, cut to 45µM, and free-floating sections were stained with above antibodies.

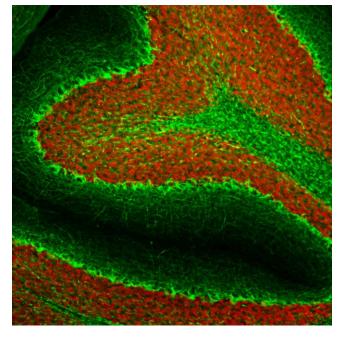

### FOR RESEARCH USE ONLY

NEUROMICS' REAGENTS ARE FOR IN VITRO AND CERTAIN NON-HUMAN IN VIVO EXPERIMENTAL USE ONLY AND NOT INTENDED FOR USE IN ANY HUMAN CLINICAL INVESTIGATION, DIAGNOSIS, PROGNOSIS, OR TREATMENT. THE ABOVE ANALYSES ARE MERELY TYPICAL GUIDES. THEY ARE NOT TO BE CONSTRUED AS BEING SPECIFICATIONS. ALL OF THE ABOVE INFORMATION IS, TO THE BEST OF OUR KNOWLEDGE, TRUE AND ACCURATE. HOWEVER, SINCE THE CONDITIONS OF USE ARE BEYOND OUR CONTROL, ALL RECOMMENDATIONS OR SUGGESTIONS ARE MADE WITHOUT GUARANTEE, EXPRESS OR IMPLIED, ON OUR PART. WE DISCLAIM ALL LIABILITY IN CONNECTION WITH THE USE OF THE INFORMATION CONTAINED HEREIN OR OTHERWISE, AND ALL SUCH RSKS ARE ASSUMED BY THE USER. WE FURTHER EXPRESSLY DISCLAIM ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE.-V2/08/2012