



## RFP Expressing Human Gastric Carcinoma N87 Cells

---

**Catalog #:** TR02-RFP

**Cell #:** >5x10<sup>5</sup> cells



**Storage:** Liquid Nitrogen until ready for culture.  
While Culturing keep in 37°C CO<sub>2</sub> incubator

**Product Format:** Frozen Vial

---

### GENERAL INFORMATION

Human Gastric Carcinoma N87 Cells (N87 cells) were established from a gastric carcinoma patient. RFP-Tagged Human Gastric Carcinoma N87 Cells (RFP-N87 cells) are selected from N87 cells after infected with lentiviruses expressing RFP with puromycin. The cells can be cultivated with Universal Full Growth Medium (TM001) and they have a minimum population doubling capacity > 12 when cultured following the detailed protocol described below.

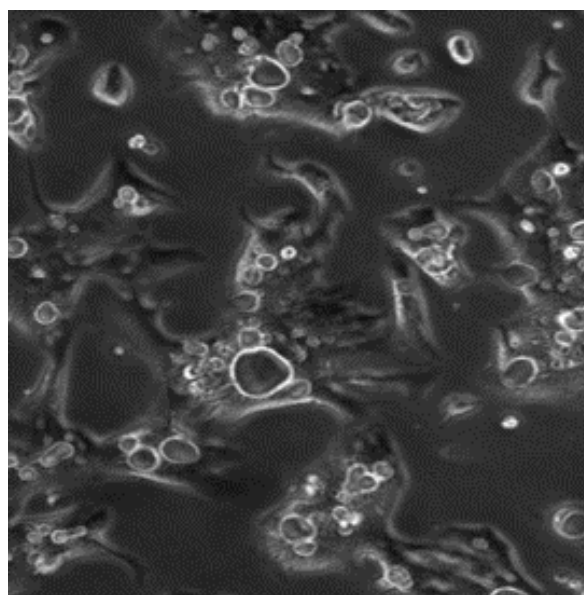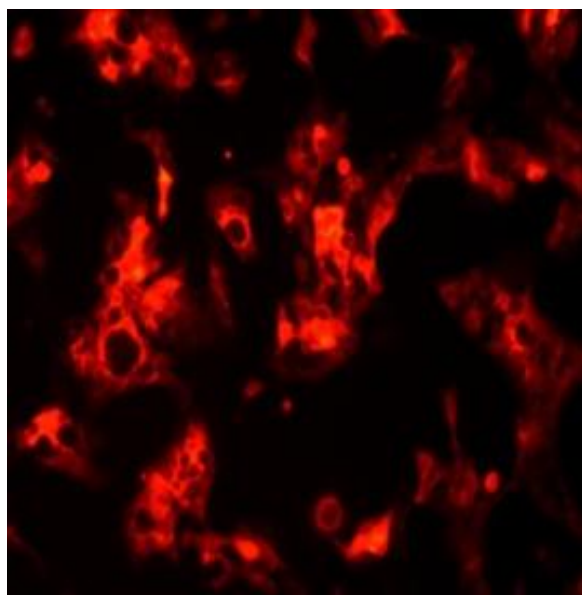

*Product is for Research use only.*

Frozen Vials are shipped in a Dry Ice Package.

### CHARACTERIZATION OF THE CELLS

N87 cells are negative for HIV-1, HBV, HCV, and mycoplasma.

### HANDLING OF ARRIVING CELLS

When you receive the dry ice package with cells in frozen vials, transfer the frozen vials of cells into a -80°C freezer for short period storage or a liquid nitrogen tank for long-term storage.

### FOR RESEARCH USE ONLY

NEUROMICS REAGENTS ARE FOR IN VITRO AND CERTAIN NON-HUMAN IN VIVO EXPERIMENTAL USE ONLY AND NOT INTENDED FOR USE IN ANY HUMAN CLINICAL INVESTIGATION, DIAGNOSIS, PROGNOSIS, OR TREATMENT. THE ABOVE ANALYSES ARE MERELY TYPICAL GUIDES. THEY ARE NOT TO BE CONSTRUED AS BEING SPECIFICATIONS. ALL OF THE ABOVE INFORMATION IS, TO THE BEST OF OUR KNOWLEDGE, TRUE AND ACCURATE. HOWEVER, SINCE THE CONDITIONS OF USE ARE BEYOND OUR CONTROL, ALL RECOMMENDATIONS OR SUGGESTIONS ARE MADE WITHOUT GUARANTEE, EXPRESS OR IMPLIED, ON OUR PART. WE DISCLAIM ALL LIABILITY IN CONNECTION WITH THE USE OF THE INFORMATION CONTAINED HEREIN OR OTHERWISE, AND ALL SUCH RISKS ARE ASSUMED BY THE USER. WE FURTHER EXPRESSLY DISCLAIM ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. v1-09809

**[www.neuromics.com](http://www.neuromics.com)**

Neuromics Antibodies • 5325 West 74<sup>th</sup> Street, Suite 8 • Edina, MN 55439  
phone 866-350-1500 • fax 612-677-3976 • e-mail [pshuster@neuromics.com](mailto:pshuster@neuromics.com)

## PROTOCOL FOR THAWING THE CELLS AND SUBCULTURE

1. Pre-coating of T25 flasks: Add 1 ml of Fibronectin (FB001) into one T25 flask and make sure whole surface of the flask is covered with the coating solution and leave flask at room temperature for overnight or 37°C for 1 hour. Dispose excessive fibronectin by aspiration. Rinsing the flask with PBS once before the flask is ready to be used.
2. Thaw the frozen cell vial in a 37°C water bath first, and then transfer the cells into the pre-coated T25 flask with 10 ml of Universal Full Growth Medium (TM001), cells usually become confluent with 5-7 days.
3. To passage the cells, rinse the cells in a T25 flask with 5 ml HBSS (Room Temperature) twice; then add 2 ml Universal Detachment Solution (RT) (AD002) into one T25 flask; gently dispose the excessive Universal Detachment Solution within 20 seconds by aspiration.
4. Leave the T25 flask with the cells at RT or 37°C for 1 min (most cells usually will detach from the surface within 1-2 mins) or monitor the cells under a microscope until most of cells become rounded up, and then gently tap the flask against the bench surface, and the cells will move on the surface of the flask when monitoring under microscope.
5. Add 5 ml of Universal Neutralization Buffer (NB001) and spin down the cells with an 800 g centrifugation for 5 mins.
6. Re-suspend the cell pellet with 10 ml or 15 ml of Universal Full Growth Medium and transfer 5 ml each into 2 or 3 pre-coated T25 flasks (for 1/2 to 1/3 subculture ratio).
7. Change the medium every 2 or 3 days and the cells usually become confluent within 7 days (when split at a 1/3 ratio).

## FOR RESEARCH USE ONLY

NEUROMICS REAGENTS ARE FOR IN VITRO AND CERTAIN NON-HUMAN IN VIVO EXPERIMENTAL USE ONLY AND NOT INTENDED FOR USE IN ANY HUMAN CLINICAL INVESTIGATION, DIAGNOSIS, PROGNOSIS, OR TREATMENT. THE ABOVE ANALYSES ARE MERELY TYPICAL GUIDES. THEY ARE NOT TO BE CONSTRUED AS BEING SPECIFICATIONS. ALL OF THE ABOVE INFORMATION IS, TO THE BEST OF OUR KNOWLEDGE, TRUE AND ACCURATE. HOWEVER, SINCE THE CONDITIONS OF USE ARE BEYOND OUR CONTROL, ALL RECOMMENDATIONS OR SUGGESTIONS ARE MADE WITHOUT GUARANTEE, EXPRESS OR IMPLIED, ON OUR PART. WE DISCLAIM ALL LIABILITY IN CONNECTION WITH THE USE OF THE INFORMATION CONTAINED HEREIN OR OTHERWISE, AND ALL SUCH RISKS ARE ASSUMED BY THE USER. WE FURTHER EXPRESSLY DISCLAIM ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. v1-09809

**[www.neuromics.com](http://www.neuromics.com)**

Neuromics Antibodies • 5325 West 74<sup>th</sup> Street, Suite 8 • Edina, MN 55439  
phone 866-350-1500 • fax 612-677-3976 • e-mail [pschuster@neuromics.com](mailto:pschuster@neuromics.com)